



Model: OM-700

USER GUIDE



PREFACE

Welcome! Now, sit back and enjoy the wonderful OM-700 massage chair. Please carefully read this user guide for the correct usage method, and pay extra attention to the caution and remarks contain herein.



Please keep this user guide handy for future reference.

CONTENTS


| | |
|-------------------------------------|----|
| Important Safety Instructions | 1 |
| □Applicable Surroundings | 1 |
| □Unsuitable Users | 1 |
| □Terms for Using | 2 |
| Functions and Features | 3 |
| Components Diagrams..... | 4 |
| □External Diagrams | 4 |
| □Internal Diagrams | 5 |
| □Outer Control Panel | 6 |
| □ Inner Control Panel | 7 |
| □LCD | 8 |
| Preparations Before Use..... | 9 |
| □Location of Chair | 9 |
| □Floor Protection | 9 |
| □Environment | 9 |
| □Grounding Connection | 9 |
| □Relocation | 9 |
| □Installation..... | 10 |
| Operating Procedures..... | 11 |
| □Power On | 11 |
| □Time Setting | 11 |
| □Chair-Back Recline | 11 |
| □Raising the Chair-Back | 11 |
| □MP3 Music Player | 12 |
| □End Massage Session | 12 |
| After Use | 13 |
| Troubleshooting..... | 14 |
| Care and Maintenance | 15 |
| □Maintenance | 15 |
| □Storing Method | 15 |
| □Inspection | 15 |
| Specifications | 16 |

IMPORTANT SAFETY INSTRUCTIONS


- Before use, please carefully read this Important Safety Instructions.
- Please follow instructions carefully to avoid any possible personal injury or damage to the product
- For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manner and should be strictly noted.

| | |
|--|---|
|  WARNING | Once operating improperly, the user may be seriously hurt. |
|  CAUTION | Once operating improperly, the user may be hurt or the product damaged. |


APPLICABLE SURROUNDINGS

| | |
|--|---|
|  CAUTION | To avoid any electrical shock or injury, do not use it in a bathroom or damp place |
| | To avoid fire, do not use product on or near any heating device or under direct sunlight. |
| | To avoid fire, do not use product on electro-thermal carpet or other heaters |
| | Make sure to place the product on level and smooth surface at all time |

UNSUITABLE USERS People suffering from the following conditions should consult with their doctor before use





| | |
|--|--|
|  PROHIBITIONS | People under treatment, injuries or abnormalities |
| | Serious tumour patients |
| | Pregnant or menstruating women |
| | Acute disease |
| | Heart attack |
| | Cardiac problems |
| | Use of heart pacemakers or any other internal medical electrical apparatus |

INSTRUCTIONS ON USING REMOTE CONTROL

| | |
|---|---|
|  CAUTION | Do not operate the remote control with wet hands |
| | Do not throw onto the ground or under water |
| | Do not pull connecting cord |
| | Do not place anything heavy on remote control |
| | It is advisable to keep the remote control on its stand |

IMPORTANT SAFETY INSTRUCTIONS

TERMS FOR USING


| | |
|---|---|
|  | When descending the footrest or backrest, use moderate pressure and be sure there are no obstructions or small children near or under the footrest. |
| | Before use, uncover back cushion and make sure inside cushion is in good order as with all other parts. Should any part be damaged contact distributor immediately and do not use until repaired. |
| | To avoid any fire or damage to the product, rated voltage must be used only. |
| | If the supply cord is damaged, the manufacturer or its service agent or similarly qualified persons must replace it in order to avoid hazard. |
| | Do not allow children or unconscious people to use product without any assistance. |
| | Do not allow children to play on the product. |
| | To avoid any injuries, do not use it immediately after meal. |
|  | Do not have more than one person on the product at any one time so as to avoid personal injury, or damaging the product. |
| | Do not use pointed or sharp objects on or near product at any time and avoid dragging on material of product. |
| | Do not use it while smoking to avoid fire. |
| | Do not use other medical apparatus when operating the product. |
| | Should any abnormalities occur when using, stop to consult with doctor. |
| | When disconnecting plug from power socket, be sure to hold insulating bush instead of power cord. |
| | Do not move product on the fragile wood or carpet floor in case any damage to the floor or truckle. |
| | Do not use it for medical purposes. |
| | Do not use petrol, caustic agents, or pesticide to clean the product. |
|  | Do not sit on backrest, armrest and footrest to avoid any injury or damage. |
| | Allow only authorized technicians to handle any breakdown, repair or alteration to the product. |
|  | When storing or maintenance, remove plug from power socket and make sure not to operate by wet hands to avoid any hazard. |
| | After use, turn power off. |

FUNCTIONS AND FEATURES

1. The stimulated chair-back mechanism is designed to automatically search and memorize the body physique and weight, and offer the best massage efficiency accordingly.



2. The chair-back mechanisms are designed to simulate the actions of an actual masseuse, and provide delicate kneading, tapping, shiatsu, and synchronization massage modes. Furthermore, the adjustable massage intensity and massage roller width help to relieve muscular fatigue and increase blood circulation.

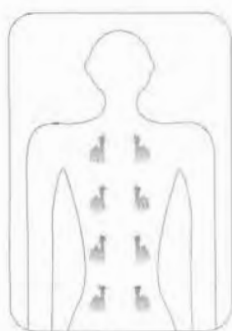


Knead



Tapping




Synchronization

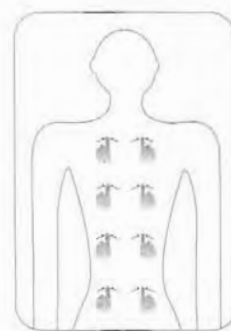


Pressing

3. The chair's seat cushion and footrest contain air bags in areas of the buttock, thighs, calves, and feet. The air bags continuously inflate and deflate to offer soothing pressure massage, and assist in relaxation of the nerves and improve blood circulation in the lower half of the body.

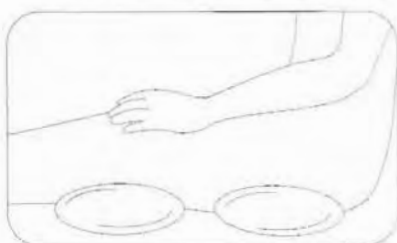

Thigh and Buttock Airbags

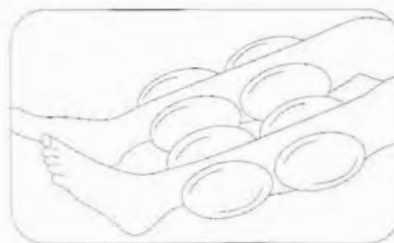

Legs and Feet Airbags

4. Automatic Overall Programs



Stops all massage functions and reposition the chair-back and footrest automatically.



Full Back AUTO: Select upper body automatic overall program.



Neck & Shoulders AUTO: Select neck & shoulders automatic partial program.



Back & Waist AUTO: Select back & waist automatic partial program.




Air Massage AUTO: Activate lower body automatic air pressure program.



Air Massage Strength: Select air pressure strength from "weak" and "strong".

COMPONENTS DIAGRAMS

EXTERNAL DIAGRAMS

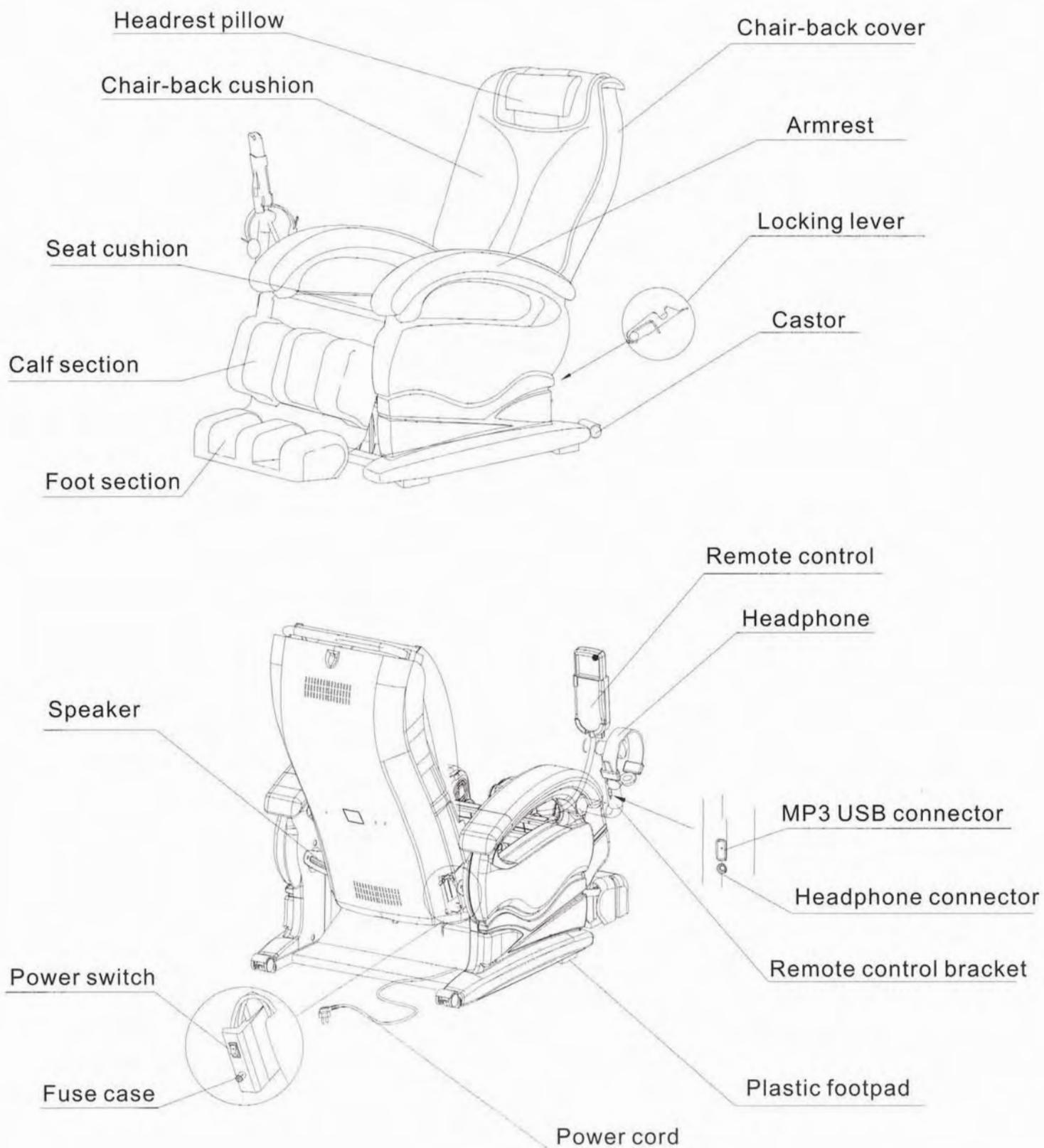

COMPONENTS DIAGRAMS


INTERNAL DIAGRAM

Chair-back massage mechanism

Seat cushion airbags




Calf airbags

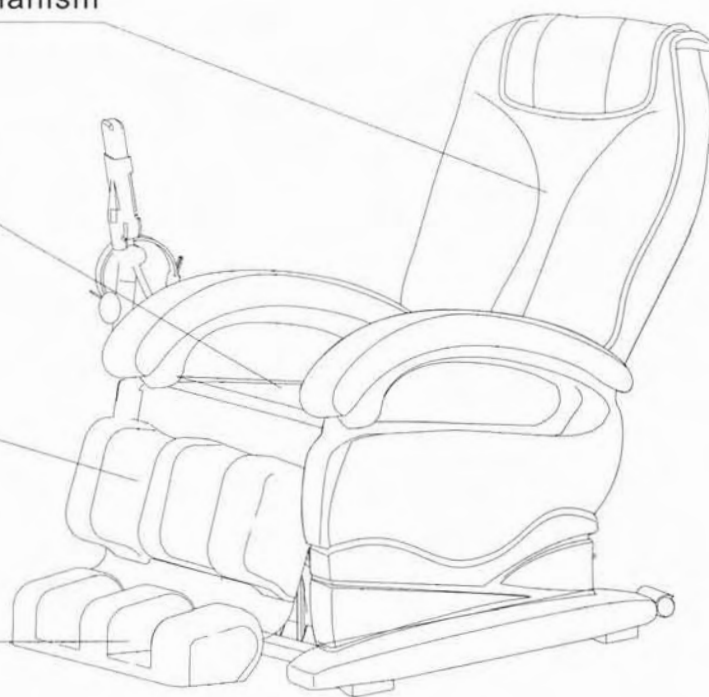Foot airbags






COMPONENTS DIAGRAMS



OUTER CONTROL PANEL

-  Turns on/off massage functions.
-  Selects massage sessions from 10, 20, or 30 minutes.
-  Stops all massage functions and reposition the chair-back and footrest automatically.





● Upper body automatic programs

-  **A** Full Back AUTO: Select upper body automatic overall program.
-  **B** Neck & Shoulders AUTO: Select neck & shoulders automatic partial program.
-  **C** Back & Waist AUTO: Select back & waist automatic partial program.

● Lower body automatic programs

-  **D** Air Massage AUTO: Activate lower body automatic air pressure program.
-  **E** Air Massage Strength: Select air pressure strength from “weak” and “strong”.








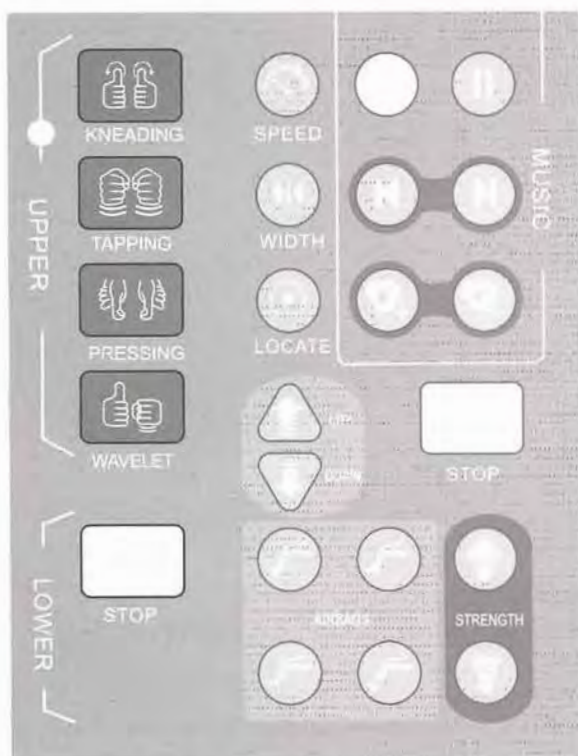
-  **Chair-back up:** Press and hold down this button to raise the chair-back and lower the footrest simultaneously. When the chair-back is completely upright, a “beep” sound is heard.
-  **Chair-back down:** Press and hold down this button to lower the chair-back and raise the footrest simultaneously. When the chair-back is completely reclined, a “beep” sound is heard.
-  **Footrest up:** Press and hold down this button to raise the footrest. When the footrest is completely raised, a “beep” sound is heard.
-  **Footrest down:** Press and hold down this button to lower the footrest. When the footrest is completely lowered, a “beep” sound is heard.

COMPONENTS DIAGRAMS






INNER CONTROL PANEL

MP3 Music Control



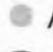







-  Play: plays music.
-  Pause: pauses music.
-  Previous: selects the previous song on the play list.
-  Next: selects the next song on the play list.
-  Volume +: increases volume of music.
-  Volume -: decreases volume of music.



Lower Body Massage


-  Air pressure massage manual control
-  Selects feet and calves air pressure massage.
-  Selects feet and seat cushion air pressure massage.
-  Selects calves and seat air pressure massage.
-  Selects feet, calves, and seat air pressure massage.

Upper Body Massage



-  Kneading: selects only the kneading massage function.
-  Tapping: selects only the tapping massage function.
-  Pressing: selects only the Shiatsu massage function.
-  Wavelet: selects both the kneading and tapping synchronize functions.
-  Speed: Selects from three chair-back massage intensities from MILD, MEDIUM, or STRONG. This adjustment is invalid when Pressing massage function is selected.
-  Width: Selects from three massage-roller widths from NARROW, REGULAR, or WIDE. This adjustment is invalid when either Kneading or Wavelet massage functions is selected.
-  Locate: Selects massage programs from Overall, Partial (Shoulders, Back, Waist), or Point.
-  Message-rollers up: Controls the massage-rollers up during Point mode.
-  Message-rollers down: Controls the massage-rollers down during Point mode.
-  Stop: Stops the chair-back massage and returns the massage-rollers to the top.

COMPONENTS DIAGRAMS

LCD

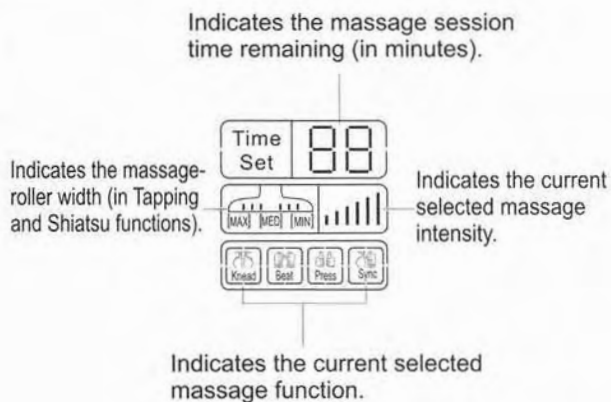


1

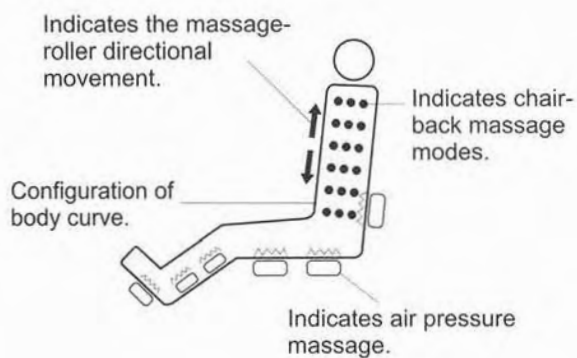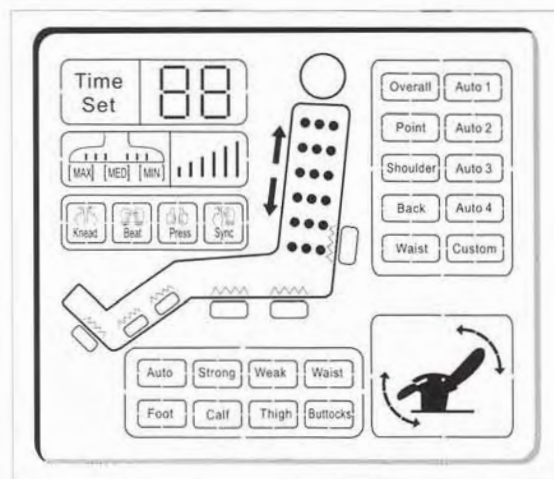


2

Indicates chair-back massage mode.

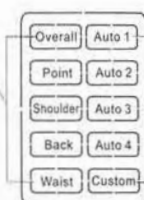

Indicates chair-back massage program.

3

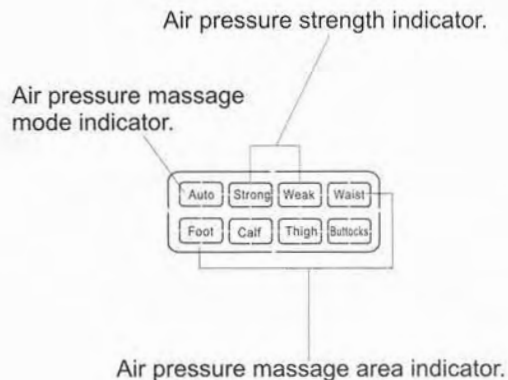

4

Chair-back up/down indicator.



Footrest up/down indicator.

5


PREPARATIONS BEFORE USE

LOCATION OF CHAIR

Safety distance around massage chair.

The massage chair-back should be at least 20 inches away from any obstacle.

The massage footrest should be at least 6 inches away from any obstacle.

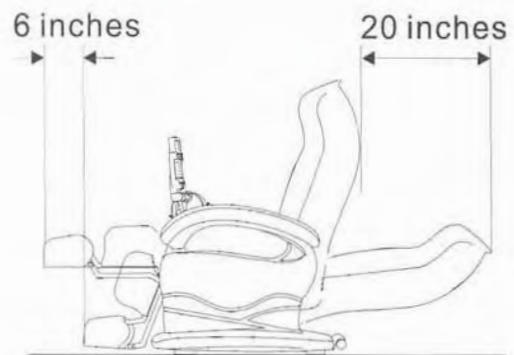

FLOOR PROTECTION

This product may cause damage to the floor, so a material (carpet for instance) under the chair to act as a medium is strongly advised.

Carpet
Floor



ENVIRONMENT



Do not keep the product in a bathroom, or a damp place, to avoid personal injury and damaging the product.



Do not keep the product near a stove, or other heating devices, to avoid fire and damage to the leather upholstery.


GROUNDING CONNECTION



This product is a classification 1 electric appliance. Please be sure to confirm ground connection to avoid electric short or accidental damage.


RELOCATION

Lift the front of the chair by the footrest, not by the seat cushion or the chair-back, and push with both hands to start moving. Do not attempt to relocate the chair while it is occupied. Do not drag or roll the chair on delicate flooring, but have two (2) or more people lift the chair completely off the ground instead. Please handle with care as the chair is very heavy.

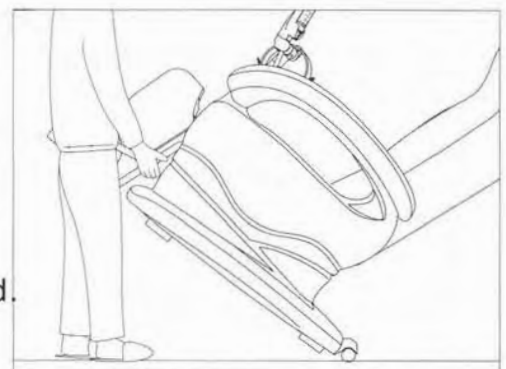


PREPARATIONS BEFORE USE

INSTALLATION

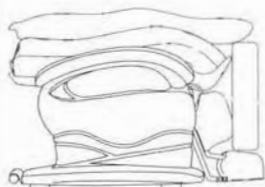

Step 1:
Hold and lift the massage chair out of its box by the points indicated in the figure. Make sure to have at least four (4) people for this step in order to avoid personal injury.


Step 2:
Situate the chair-back to its upright position until a “click” sound is heard indicating it is locked in place.

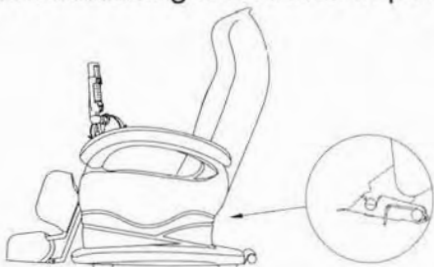


Step 3:

Connect the color-coded air tube and Velcro the seat cushion onto the frame of the chair.

Windpipe

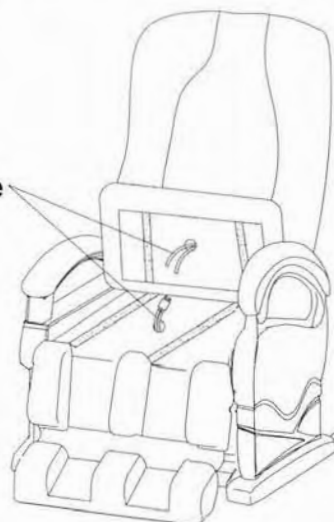


Step 4:
Velcro the headrest pillows onto the chair-back.

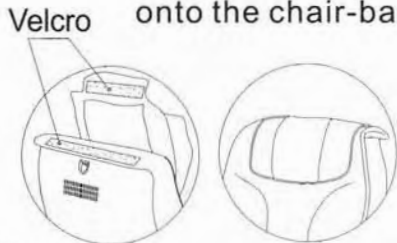

Step 5:

Remove the four (4) screws on the right armrest where the remote control holder will be attached. Slide the excess wires into the hole on the right armrest and secure the remote control holder onto the right armrest with the four (4) screws.

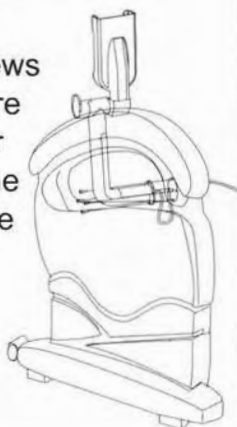

Step 6:

Connect the power cord into the wall outlet, and turn the main power switch located behind the right armrest on. The chair is ready for use.



OPERATING PROCEDURES



POWER ON


- Connect the chair's power cord to the wall outlet.
- Turn the chair's main power switch (located behind the right armrest) on.

TIME SETTING


Choose massage sessions of 10, 20, or 30 minutes. When the timer ends, the massage session will terminate. Press the Settle button to return the chair to its upright position. The default time of a massage session is 20 minutes.

CHAIR-BACK RECLINE

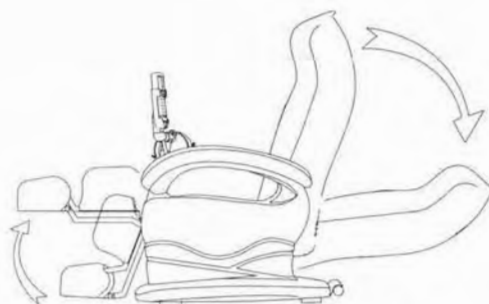

Press and hold the  button to recline the chair-back. This automatically raises the footrest simultaneously. Release the  button to stop the lowering and rising the chair-back and footrest, respectively. If a beeping sound is heard, the chair-back and/or the footrest has reached its limit.


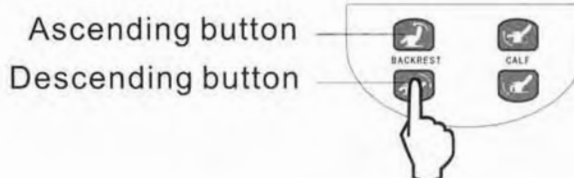


Reclining the chair-back



RAISING THE CHAIR-BACK

Press and hold the  button to raise the chair-back. This automatically lowers the footrest simultaneously. Release the  button to stop the rising and lowering of the chair-back and footrest, respectively. If a beeping sound is heard, the chair-back and/or the footrest has reached its limits.

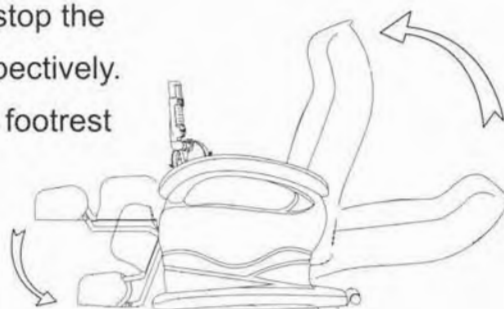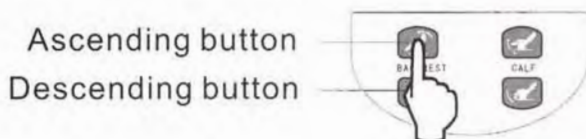

Raising the chair-back

OPERATING PROCEDURES

MP3 MUSIC PLAYER

This product includes a built-in MP3 music player.

● Storing and deleting MP3 music


❗ Please shut off the chair's main power switch (located behind the right armrest) before performing the following:

1. Connect the MP3 music player with a computer having Windows 2000 or better, using an USB cable.
2. Computer will automatically search and detect this plug and play device.
3. Open up My Computer and search for the new icon MP3 drive.
4. You can then transfer and delete songs to and from the MP3 drive.
5. Unplug the USB cable from the computer and the chair before turning the chair's power on.

For instructions on how to operate the MP3 music player, please refer to the MP3 Music Control section of Components Diagrams on page 7.

END MESSAGE SESSION

- The massage session ends when the timer reaches zero (0), or when the power button on the remote control is pressed.
- Press the settle button on the remote control will also achieve the same result as pressing the power button.
- The chair-back and footrest will only return to its upright position when the power or settle buttons on the remote control is pressed.

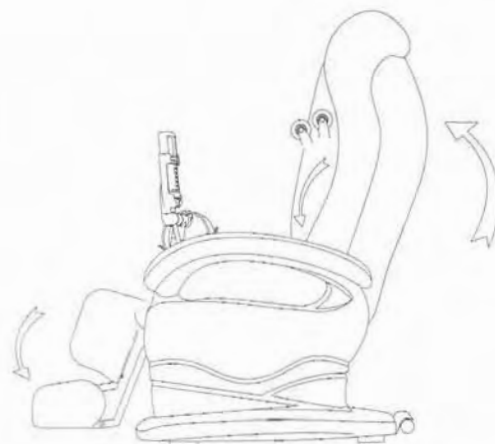


AFTER USE...

- After use, it is suggested that the chair's power switch should be turned off, and the power cord unplugged from the wall outlet.
- Keep the remote control in its holder.

Should the product not be used for a long period of time, fold the chair-back by pushing down on the lock lever located beneath the chair-back (Refer to page 4), and store the product in a cool and dry environment.

Storage

Turn the chair's main power switch off, unplug its power cord from the wall outlet, and fold the chair-back downwards by pushing down on the lock lever located beneath the chair-back. Store the chair in a cool and dry place.

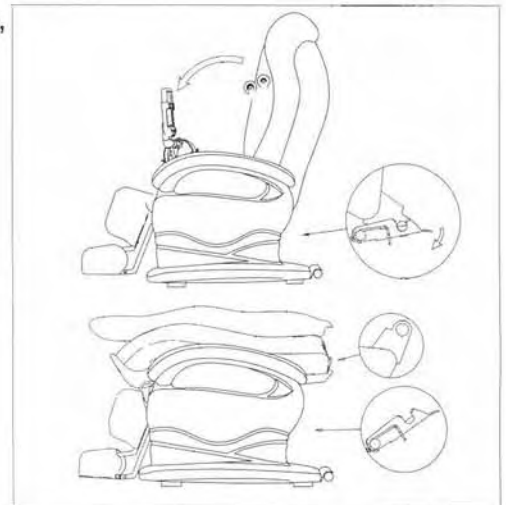

TROUBLESHOOTING

| Symptoms | Please check | Solutions |
|---|---|--|
| Noise from the chair-back during massage session. | None | Normal working condition. |
| Noise from air compressor. | None | Normal working condition. |
| Not operating. | If the fuse is blown. | Replace fuse of same model. |
| | If the chair's main power switch is turned off. | Turn the chair's main power switch on. |
| | If power cord became disconnected from the wall outlet. | Reconnect power cord to the wall outlet. |
| Abnormal playing of certain music. | If music is in MP3 format. | Replace music format. |
| Computer cannot recognize MP3 music player. | If the chair's main power switch is turned off. | Turn the chair's main power switch on and try again. |
| | If PC uses Windows 2000 or better operating system. | Use a PC that has Windows 2000 or better operating system. |
| No volume from the headphone. | If the headphone is connected to the remote control holder. | Unplug and reinsert the headphone wire. |
| | If the volume is turned down or off. | Turn the volume up. |
| Motors overload. | If too much pressure is placed against the chair-back and footrest. | Decrease the amount of pressure against the chair-back and footrest. |

If the problem persists, please turn off the chair, unplug its power cord from the wall outlet, and contact the customer support for guidance.

CARE AND MAINTENANCE

Caution Always unplug the chair from the wall outlet before performing maintenance on the chair. Do not attempt to unplug the chair with a damp cloth or wet hands.

MAINTENANCE

Main body

To remove dirt from plastic parts and armrests, wipe with a damp cloth using mild cleaning agent and dry thoroughly with a clean and dry cloth.

NOTICE



Do not use alcohol, caustic agents, detergents, and etc. to clean the chair so as to avoid damage, discoloration, cracking, and etc. Of the upholstery.

Remote control

Wipe remote control with a dry cloth only.

NOTICE



Never use damp materials to clean remote control to avoid accidents.

Headrest pillow and chair-back cushion

If dirty, clean with a damp cloth and mild cleaning agent, followed with a clean damp cloth, and air dry.

STORING METHOD

After dusting and cleaning, keep the chair in a cool and dry place. If storing for a long period of time, cover the chair.

NOTICE



Avoid placing the chair in direct sunlight, or near excessive heat, as this will cause fading and discoloration.

INSPECTION



Please inspect, and stop using if...

- There is a severe thunderstorm,
- There is a parched smell,
- You suspect an electrical wire overheating, or
- There are other abnormal phenomena.

Stop using

In case of accident, please pull off the plug and contact agent or our company. Do not allowed to disassemble or fix by yourself.

SPECIFICATIONS

Model: OM-700

Voltage: 120V~ 60Hz

Power:175W

Dimensions:35"x 51" x 47" (W x D x H) - Upright

35" x 70" x 33" (W x D x H) - Fully Reclined

Net Weight: 188 lbs.

Standardization:Class 1

ALL RIGHTS RESERVED, ANY DISCREPANCIES PLEASE REFER TO
THE ACTUAL PRODUCT.



OMEGA MASSAGE CHAIR, INC.

Tel: 1-800-659-3650

www.omegamassage.com

Email Address: info@omegamassage.com